

Catheter Production Fixtures

...and other extras



1 • Balloon Trimmer/Tip Cutter

Performs precision square end cuts to specified lengths. Quick and repeatable results on items such as balloon necks and catheter tips.




2 • Slitting Fixture, Pneumatic or Manual

This fixture is used to make axial slits precisely aligned with the tubing centerline through one wall only - such as for forming an inner lumen opening on PTCA distal shaft tubing. The pneumatic model features a foot pedal control.



3 • Marker Band Pre-Crimper

The Marker Band Pre-Crimper simultaneously positions and pre-crimps marker bands on a catheter shaft in preparation for the swaging process.

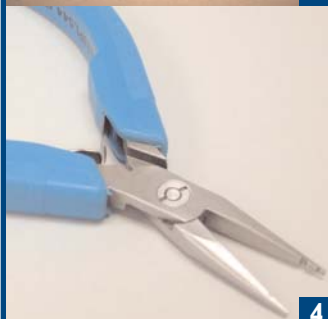

4 • Crimping Pliers

The Crimping Pliers are used to manually position marker bands on the catheter shaft. The pliers are custom made to your specifications, with up to 5 crimping grooves depending on the sizes required.

5 • Catheter Stand

This rollaway stand is ideal for transferring semi-finished and finished catheters between work stations. The stand holds approximately 80 catheters. Catheters are hanging by the luer hub that is retained in a slot on the side of the cross bars. Placement and removal of catheters is fast and easy.



27752 El Lazo Road • Laguna Niguel, CA 92677 • Phone: 949.448.7056 • Fax: 949.448.7016
www.interfaceusa.com • info@interfaceusa.com